



airMmax

Perfection has a name



airMax

Air polishing by Satelec®

*Fast and straightforward prophylactic treatment for perfect results every time
With AirMax, choose simplicity, performance, comfort and efficiency*



Built-in water heater

Thanks to its water heating system, the water/ powder mix is at a comfortable temperature for the patient

Automatic purge

One press on the Purge button will flush out the device's complete internal system and avoid the risk of untimely clogging

Transparent powder container

The transparent container makes filling easier and the remaining powder level can be checked at a glance



Plug and Play

Versatile and mobile

The simple quick-coupling connection allows leak-free installation of AirMax, just like a mobile unit, in different treatment rooms

Prophy Spray

- No more clogged handpieces...
- No more returns to After Sales for repairs...
- ProphySpray is the first ever maintenance agent for air polishers

Fine, regular and thorough polishing for a sparkling smile!

Innovative design means an end to clogging


The special design of the ProphyPen handpiece and its titanium nozzle keeps the water completely separate from the air-powder jet right up to the last second. Clogging is now a thing of the past!

Easy maintenance means faster treatment

The nozzle and the ProphyPen handpiece can be disconnected - it only takes a few seconds to replace soiled elements before treating the next patient

Priority to asepsis means more safety

The ProphyPen is entirely autoclavable to keep the risk of cross-contamination between patients at bay

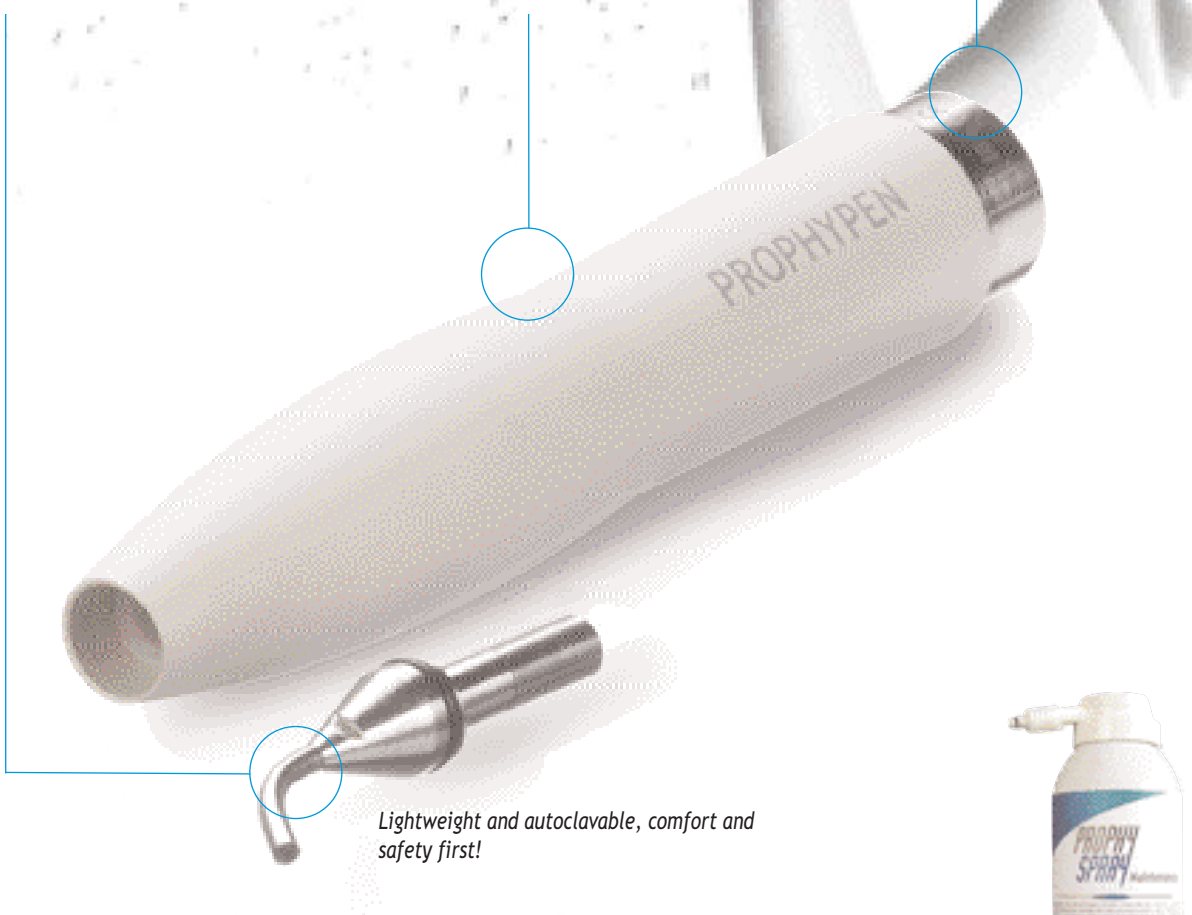

Lightweight and autoclavable, comfort and safety first!

ProphySpray dissolves any remaining traces of powder and guarantees the smooth flow and maximum efficiency of the water/air/powder mix:

- The ProphyPen nozzle maintains its efficiency for longer

- It takes effect immediately, with no hitches and no risks: made from only 100% natural constituents
- ProphySpray presents no danger to your patients or practice personnel and cannot cause any material damage to your device
- Treatments are carried out as they should be: smoothly, and giving a better return on prophylactics



airMax

- Fine, regular and thorough polishing
- Hardwearing and reliable
- Greatly reduced risk of clogging
- Work in greater comfort, faster and without hand fatigue
- “Plug and Play” quick-coupling installation



Powder Max

Kind on enamel

The powder's formula ensures perfect polishing and respect of the enamel

No risk of clogging

Its grit size has been selected to guarantee efficiency but avoid all risk of clogging

Just use what you need

Measuring the quantity of powder needed to fill the container is easy thanks to the bottle cap

Available in spearmint or lemon flavour

Technical characteristics

Power supply: 115 VAC / 230 VAC
Water pressure: 1 to 5 Bar
Air pressure: 4.5 to 6 Bar

Dimensions:
Width 250 mm
Height 110 mm
Depth 300 mm

The nozzle and handpiece are sterilizable in a class B autoclave according to EN 13060 Standard:
Temperature: 134°C ;
Pressure : 2 Bars (200 hPa) ;
Sterilization time: 18 minutes under steam heat
See individual product instructions manuals for further information

This medical device is referenced as class IIa according to the current European Directive in force
This medical device is manufactured in compliance with the current Standard (IEC 60601-1)
This medical device has been designed and is manufactured according to the EN ISO 13485 quality assurance system requirements.
Please refer to specific Standards and regulations in force in your country

